


GOPHER™ **GOLD** Catheter

*A new twist
on crossing stenoses*



- Tapered threaded distal tip engages stenosis to support passage
- Gold-plated tip for bright radiopacity
- Tightly braided backbone for enhanced torque without breakage
- Sealed polymer jacket allows infusion of medication and contrast
- Pure PTFE lining to facilitate wire movement

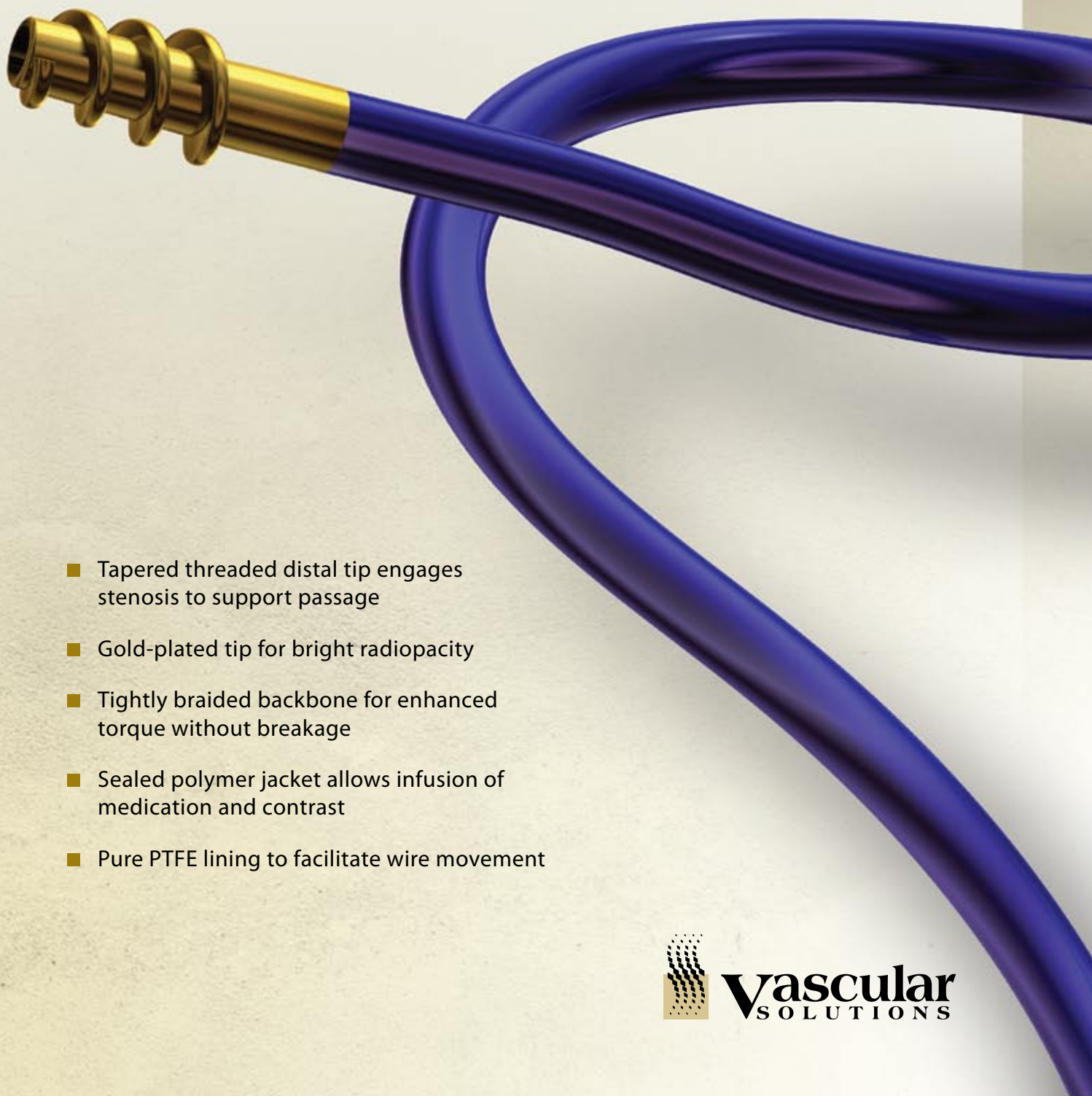

vascular
SOLUTIONS

GOPHERTM **GOLD** Catheter



Clockwise rotation of the Gopher Gold catheter while in contact with the stenosis anchors the tip and increases support.



Further clockwise rotation *pulls* the Gopher Gold catheter through the stenosis.




The Gopher Gold catheter's unique tapered threaded distal tip design provides rotational support for crossing stenoses.

Tapered, threaded distal tip translates clockwise rotation into support for passage

- Tapered tip for engaging stenosis
- Clockwise rotation pulls catheter through stenosis, counter-clockwise withdraws it
- Smoothly polished, non-augering threads
- Highly visible gold radiopaque tip

Multilayer construction provides both flexibility and torque transmission

- Proximal nitinol sleeve transfers 1:1 torque during rotation
- Braided shaft provides outstanding flexibility for wire tracking
- PTFE liner provides smooth wire delivery
- Sealed polymer outer jacket for delivering medication and contrast infusion

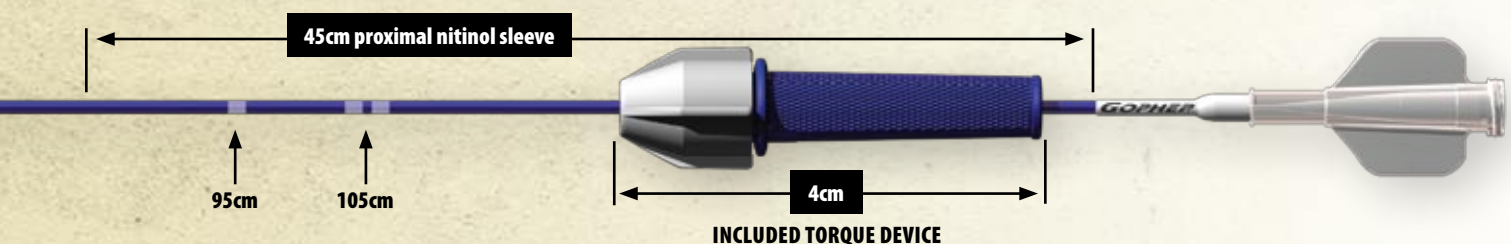

GOPHER™ GOLD Catheter


Model	Description	Box Quantity
5620	Gopher Gold Catheter	1

Gopher Gold Specifications

Guide catheter compatibility:	≥5F (0.058"/1.47mm min. I.D.)
Guidewire compatibility:	≤0.014"
Distal tip O.D.:	1.92F/0.025" (0.635mm)
Tip length:	0.115" (2.92mm)
Crossing profile:	3F/0.039" (1mm)
Working length:	135cm
Positioning marks:	95cm = single white mark 105cm = double white marks
Radiopaque marker:	Threaded tip

The Gopher Gold catheter is intended to be used in conjunction with steerable guidewires to access discreet regions of the coronary and/or peripheral vasculature. It may be used to facilitate placement and exchange of guidewires and other interventional devices and to subselectively infuse/deliver diagnostic and therapeutic agents.

The Gopher Gold catheter is contraindicated for high-pressure injections and for use in the cerebral vasculature.

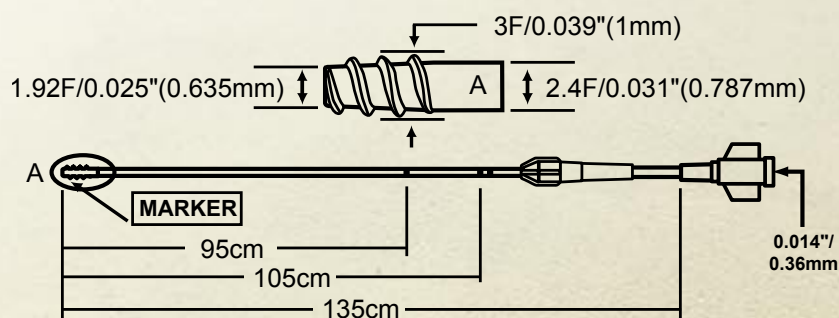

Vascular Solutions, Inc.
6464 Sycamore Court
Minneapolis, Minnesota 55369 USA

Customer Service:
United States: 888.240.6001
International: (001) 763.656.4298
customerservice@vascularsolutions.com



www.vascularsolutions.com

CAUTION: Federal Law (U.S.A.) restricts this device to sale by or on the order of a physician. Please see the Instructions for Use for a complete listing of the indications, contraindications, warnings and precautions.
U.S. and International patents pending.
Gopher is a trademark of Vascular Solutions, Inc.
©2010 Vascular Solutions, Inc. All rights reserved.
ML2105 Rev. B 01/10