


REFERENCES—all data/analysis on file.*


1. Proximal Shaft Stability – Gurley Testing (N = 5) Echelon 10, N = 5 Excelsior SL-10[™]), data/analysis on file.

Catheter type	Avg. prox. shaft stiffness (mg) at 15cm from the hub	Std. Dev. Prox. shaft stiffness (mg) 15cm from the hub	% stiffer
Echelon 10	530	29	45 %
Excelsior SL-10	365	21	control

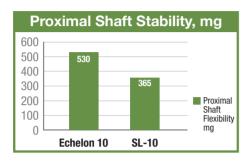

2. Distal Tip Flexibility – Tinius Olsen Testing (N = 5 Echelon 10, N = 3 Excelsior SL-10), data/analysis on file.

Catheter type	Avg. distal tip stiffness at 2mm from distal tip (10 ⁴ , inch-lbs)	St. Dev. distal tip stiffness at 2mm from distal tip (10 ⁴ in-lbs)	% flexible
Echelon 10	1.8	0.2	54 %
Excelsior SL-10	3.9	0.3	control

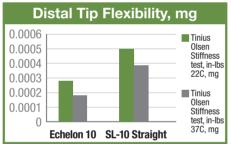

3. Friction Protection ev3 InnerGlide™ layer of advanced Pebax coating between PTFE and Progressive Nitinol Braid helps facilitate smooth, predictable responsiveness.



4. Torque Response Testing - Summary of Torque Response (N = 5Echelon 10, N = 5 Excelsion SL-10), data/analysis on fil

	Catheter type	Average # of turns for torque response	
r ile.	Echelon 10	6.2	
	Excelsior SL-10	16.4	

5. Crush Test

Laboratory testing image of Echelon 10 and Excelsior SL-10 catheters after having "alligator" clamps applied for approximately 5 minutes, then removed. Data/analysis on file.



SL-10 after crush test

Echelon 10 after crush test

Echelon and InnerGuide are trademarks of ev3. All other marks are the property of their respective owners. Indications, contraindications, warnings and instructions for use can be found in the product labeling supplied with each device. * Data on file. FD-2591

NEUROVASCULAR | PERIPHERAL VASCULAR Access · Balloons · Embolic Coils · Embolic Protection · Flow Diversion · Flow Restoration · Liquid Embolics · Plaque Excision · Procedural Support · Remodeling Devices · Retrieval Devices · Stents ev3 Neurovascular 9775 Toledo Way Irvine, CA 92618 USA PH +1 949 837 3700 ev3 S.r.l Italy PH +39 0267 977 61 FX +39 0266 711 637 ev3 GmbH Germany, Austria PH +49 228 528 830 ev3 Europe International Headquarters 106-108 rue La Boétie 75008 Paris ev3 SAS France PH +33 (0) 156 88 31 10 FX +33 (0) 156 88 31 11 FX +49 228 528 8360

Plymouth, MN 55441 USA **FX** +1 949 837 2044

The Netherlands **PH** +31 (0) 433 659 220 **FX** +31 (0) 43 364 6395

ev3 Technologies Iberica, S.L. Spain PH +34 91 656 7154 FX +34 91 656 7214


ev3 B.V. Benelux PH +31 (0) 433 659 223 FX +31 (0) 433 650 283 ev3 Nordic AB PH +46 859 000 950 FX +46 859 000 959

ev3 Sp z o.o. Poland PH +48 32 747 01 44 FX +48 32 747 01 45

ev3 Ltd. United Kingdom PH +44 1279 659 900 FX +44 1279 654 900



www.ev3.net

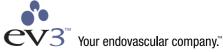

Echelon* Micro Catheters Looked into Echelon™ lately? See how much it delivers.

©2010 ev3 Inc. All rights reserved.

PH +33 156 88 59 10 FX +33 156 88 59 11

75557-002 (A) DEC 10 - Intl

Start with Echelon™. Finish with confidence.

Every procedure is different. And the tools you choose can make a difference. So we make sure Echelon micro catheters are designed to outperform the others. Need proof? Read on. Then be sure to insist on Echelon.

Echelon Micro Catheters
All catheters
are not created equal.
Here's proof.



Echelon 10 provided 45% greater proximal column stability than Excelsior SL10 measured at 15cm from hub.¹

54% MORE

The distal tip of Echelon 10 is 54% more flexible than that of the Excelsior SL10 measured at 2mm from tip.²

MORE Precision

Proprietary InnerGlide[™] layer of smooth coating between braid and PTFE helps reduce friction and increase "as-expected" performance.³

62% MORE

Echelon 10 delivered 62% better torque response in tests measuring the average number of turns required to acctuate.⁴

Visibly MORE

Post-testing images show that Echelon 10 exhibited less deformation than Excelsior SL10 after crush tests.⁵

····· Proximal Column Stability

···· Distal
Tip Flexibility

Friction

Protection

Torque Response

Crush
Protection

