

Helena Laboratories Point of Care



The MINI II is the simplest, least expensive instrument available, for activated clotting time testing, yet it retains high sensitivity, 2-point clot detection for greatest reliability. Performance, linearity and CVs parallel larger Actalyke analyzers. It's ideal for use in the catheterization or the surgical suite. Optional battery pack available.



Actalyke® MINI II ACT Analyzer

- **6" Footprint** - The MINI II features a very small footprint, but has a very large high visibility screen. Operation is extremely simple. Just insert an ACT reagent tube containing whole blood and push the button.
- **Full Cross-Compatibility** - The MINI II is fully compatible with all ACT tubes. Helena offers a full line of Actalyke ACT tubes, including traditional glass bead, kaolin and celite, plus new MAX-ACT.

Actalyke MINI II Specifications

Tests: Whole Blood Activated Clotting Time (ACT). Accepts Actalyke ACT Tubes with celite, kaolin, glass-bead activators; MAX-ACT Tubes with blended activators; all International Technidyne Hemochron® tubes

Measurement Range: 0-1500 seconds

Incubation Temperature: 36.5°C - 37.5°C (97.7°F - 99.5°F)

Dimensions: 15.5 x 12.2 x 16.0 cm (6.1"W x 4.8"H x 6.3"L)

Display: 8 x 2.54 cm (3.2" x 1")

Weight (with printer): 2.4 kg (5.3 lbs.)

Line Voltage: 47-63 Hz, 85-265 V A/C (Battery optional - DC)

Environment: 15°C to 30°C (59 to 86°F)

Data Transfer Capability: Printer option available for hard copy

CLIA Classification: Non-waived

Actalyke MINI II and Supplies

Cat. No.	Item
5753	Actalyke MINI II Analyzer, 110/220V (w/printer)
5755	Actalyke MINI II Analyzer, 110/220V (without printer)
5763	Actalyke MINI II, 110/220V (with printer and battery)
5765	Actalyke MINI II, 110/220V (with battery)
1197	Comprehensive Replacement for Actalyke MINI II Analyzer
9134	Service Manual for Actalyke MINI II
5754	Actalyke MINI II Battery Pack
5757	Actalyke Thermometer
A-PPR	Actalyke Thermal Printer Paper (3 rolls)
XL-ECT	Actalyke Electronic Clotting Tube

Helena Laboratories Point of Care





Actalyke Reagent Tubes for Activated Clotting Time

MAX-ACT™ Tubes

Cat. No. Item

MAX-ACT MAX-ACT Tubes (Gray Flip-Top), 50 tubes/box

MAX-ACT Tubes, with an upper linearity limit of 6 units of heparin per mL of blood, are ideal for monitoring patients receiving moderate to high heparin dosages. For maximum heparin specificity, MAX-ACT contains multiple particulate activators - kaolin, celite, and glass beads - that ensure maximal Factor XII activation and compensate for potential patient sensitivity to a single activator.

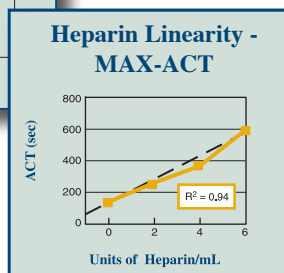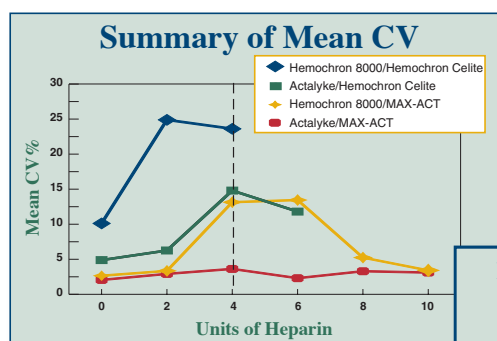

Celite Tubes

Cat. No. Item

C-ACT Celite ACT Tubes (Black Flip-Top), 50 tubes/box

C-ACT Tubes contain a celite (diatomaceous earth) activator and can be used to monitor patients undergoing procedures requiring administration of moderate to high heparin concentrations.

Kaolin ACT Tubes

Cat. No. Item

K-ACT Kaolin ACT Tubes (Yellow Flip-Top), 50 tubes/box

Kaolin ACT Tubes are useful in monitoring patients receiving moderate to high heparin concentrations in the presence of aprotinin.

Glass Bead ACT Tubes

Cat. No. Item

G-ACT Glass Bead ACT Tubes (White Flip-Top), 50 tubes/box

The Glass Bead ACT is useful for monitoring patients undergoing procedures requiring low to moderate heparin dosages.

Actalyke Quality Control

Whole Blood Control Kit for MAX-ACT, C-ACT, K-ACT

Cat. No. Item

AQC-H Whole Blood QC Kit

1 x 3 mL QC Level - I 1 x 3 mL QC Level - III
1 x 7 mL Distilled H₂O 1 x 7 mL CaCl₂ (0.02 M)

The kit provides two levels of quality control, normal and high. The lyophilized control is stable 4 hours at room temperature once reconstituted.

Whole Blood Control Kit for G-ACT

Cat. No. Item

AQC-L G-ACT Whole Blood QC Kit

1 x 3 mL QC Level - I 1 x 3 mL QC Level - II
1 x 7 mL Distilled H₂O 1 x 7 mL CaCl₂ (0.02 M)

The kit provides two levels of quality control, normal and low. The lyophilized control is stable 4 hours at room temperature once reconstituted.

Actalyke Electronic Clotting Tube

Cat. No. Item

XL-ECT Actalyke Electronic Clotting Tube

The Electronic Clotting Tube simulates patient sample clots at three therapeutic heparin levels: 100, 300 and 500 seconds for daily QC of all Actalyke instruments.