

Protégé™ EverFlex™  
Self-Expanding Peripheral Stent System  
ORDERING INFORMATION

Product Number Catheter Length 80 cm	Product Number Catheter Length 120 cm	Stent Diameter (mm)	Stent Length (mm)	Lumen Size (mm)	Sheath Size (Fr)	Recommended Guidewire (in)	Crossing Profile (in)
PRP35-05-020-080	PRP35-05-020-120	5	20	3.5 - 4.5	6	0.035	0.079
PRP35-05-030-080	PRP35-05-030-120	5	30	3.5 - 4.5	6	0.035	0.079
PRP35-05-040-080	PRP35-05-040-120	5	40	3.5 - 4.5	6	0.035	0.079
PRP35-05-060-080	PRP35-05-060-120	5	60	3.5 - 4.5	6	0.035	0.079
PRP35-05-080-080	PRP35-05-080-120	5	80	3.5 - 4.5	6	0.035	0.079
PRP35-05-100-080	PRP35-05-100-120	5	100	3.5 - 4.5	6	0.035	0.079
PRP35-05-120-080	PRP35-05-120-120	5	120	3.5 - 4.5	6	0.035	0.079
PRP35-05-150-080	PRP35-05-150-120	5	150	3.5 - 4.5	6	0.035	0.079
PRP35-06-020-080	PRP35-06-020-120	6	20	4.5 - 5.5	6	0.035	0.079
PRP35-06-030-080	PRP35-06-030-120	6	30	4.5 - 5.5	6	0.035	0.079
PRP35-06-040-080	PRP35-06-040-120	6	40	4.5 - 5.5	6	0.035	0.079
PRP35-06-060-080	PRP35-06-060-120	6	60	4.5 - 5.5	6	0.035	0.079
PRP35-06-080-080	PRP35-06-080-120	6	80	4.5 - 5.5	6	0.035	0.079
PRP35-06-100-080	PRP35-06-100-120	6	100	4.5 - 5.5	6	0.035	0.079
PRP35-06-120-080	PRP35-06-120-120	6	120	4.5 - 5.5	6	0.035	0.079
PRP35-06-150-080	PRP35-06-150-120	6	150	4.5 - 5.5	6	0.035	0.079
-	PRP35DR-06-200-120	6	200	4.5 - 5.5	6	0.035	0.079
PRP35-07-020-080	PRP35-07-020-120	7	20	5.5 - 6.5	6	0.035	0.079
PRP35-07-030-080	PRP35-07-030-120	7	30	5.5 - 6.5	6	0.035	0.079
PRP35-07-040-080	PRP35-07-040-120	7	40	5.5 - 6.5	6	0.035	0.079
PRP35-07-060-080	PRP35-07-060-120	7	60	5.5 - 6.5	6	0.035	0.079
PRP35-07-080-080	PRP35-07-080-120	7	80	5.5 - 6.5	6	0.035	0.079
PRP35-07-100-080	PRP35-07-100-120	7	100	5.5 - 6.5	6	0.035	0.079
PRP35-07-120-080	PRP35-07-120-120	7	120	5.5 - 6.5	6	0.035	0.079
PRP35-07-150-080	PRP35-07-150-120	7	150	5.5 - 6.5	6	0.035	0.079
-	PRP35DR-07-200-120	7	200	5.5 - 6.5	6	0.035	0.079
PRP35-08-020-080	PRP35-08-020-120	8	20	6.5 - 7.5	6	0.035	0.079
PRP35-08-030-080	PRP35-08-030-120	8	30	6.5 - 7.5	6	0.035	0.079
PRP35-08-040-080	PRP35-08-040-120	8	40	6.5 - 7.5	6	0.035	0.079
PRP35-08-060-080	PRP35-08-060-120	8	60	6.5 - 7.5	6	0.035	0.079
PRP35-08-080-080	PRP35-08-080-120	8	80	6.5 - 7.5	6	0.035	0.079
PRP35-08-100-080	PRP35-08-100-120	8	100	6.5 - 7.5	6	0.035	0.079
PRP35-08-120-080	PRP35-08-120-120	8	120	6.5 - 7.5	6	0.035	0.079
PRP35-08-150-080	PRP35-08-150-120	8	150	6.5 - 7.5	6	0.035	0.079
-	PRP35DR-08-200-120	8	200	6.5 - 7.5	6	0.035	0.079

INDICATIONS: The Protégé EverFlex Stent is indicated for use in common iliac, external iliac, superficial femoral, proximal popliteal, and subclavian arteries.  
Protégé EverFlex is a trademark of ev3 Inc.  
Indications, contraindications, warnings and instructions for use can be found in the product labeling supplied with each device.

NEUROVASCULAR | PERIPHERAL VASCULAR

Access · Balloons · Embolic Coils · Embolic Protection · Flow Diversion · Flow Restoration · Liquid Embolics · Plaque Excision · Procedural Support · Remodeling Devices · Retrieval Devices · Stents

ev3 Europe  
International Headquarters  
106-108 rue La Boétie  
75008 Paris  
France  
PH +33 156 88 59 10  
FX +33 156 88 59 11

ev3 Corporate  
World Headquarters  
Peripheral Vascular  
3033 Campus Drive  
Plymouth, MN 55441  
USA  
PH +1 763 398 7000  
FX +1 763 398 7001

ev3 Neurovascular  
9775 Toledo Way  
Irvine, CA 92618 USA  
PH +1 949 837 3700  
FX +1 949 837 2044  
  
ev3 International  
Distribution Centre  
Europalaan 25  
6199 AB Maastricht-Airport  
The Netherlands  
PH +31 (0) 433 659 220  
FX +31 (0) 43 364 6395

ev3 SAS France  
PH +33 (0) 156 88 31 10  
FX +33 (0) 156 88 31 11  
  
ev3 B.V. Benelux  
PH +31 (0) 433 659 223  
FX +31 (0) 433 650 283

ev3 Technologies  
Iberica, S.L. Spain  
PH +34 91 656 7154  
FX +34 91 656 7214

ev3 S.r.l. Italy  
PH +39 0267 977 61  
FX +39 0266 711 637

ev3 Nordic AB  
PH +46 859 000 950  
FX +46 859 000 959

ev3 Sp z o.o. Poland  
PH +48 32 747 01 44  
FX +48 32 747 01 45


ev3 GmbH Germany, Austria  
PH +49 228 528 830  
FX +49 228 528 8360

ev3 Ltd. United Kingdom  
PH +44 1279 659 900  
FX +44 1279 654 900

www.ev3.net



EverFlex™  
Self-Expanding  
Peripheral Stent System

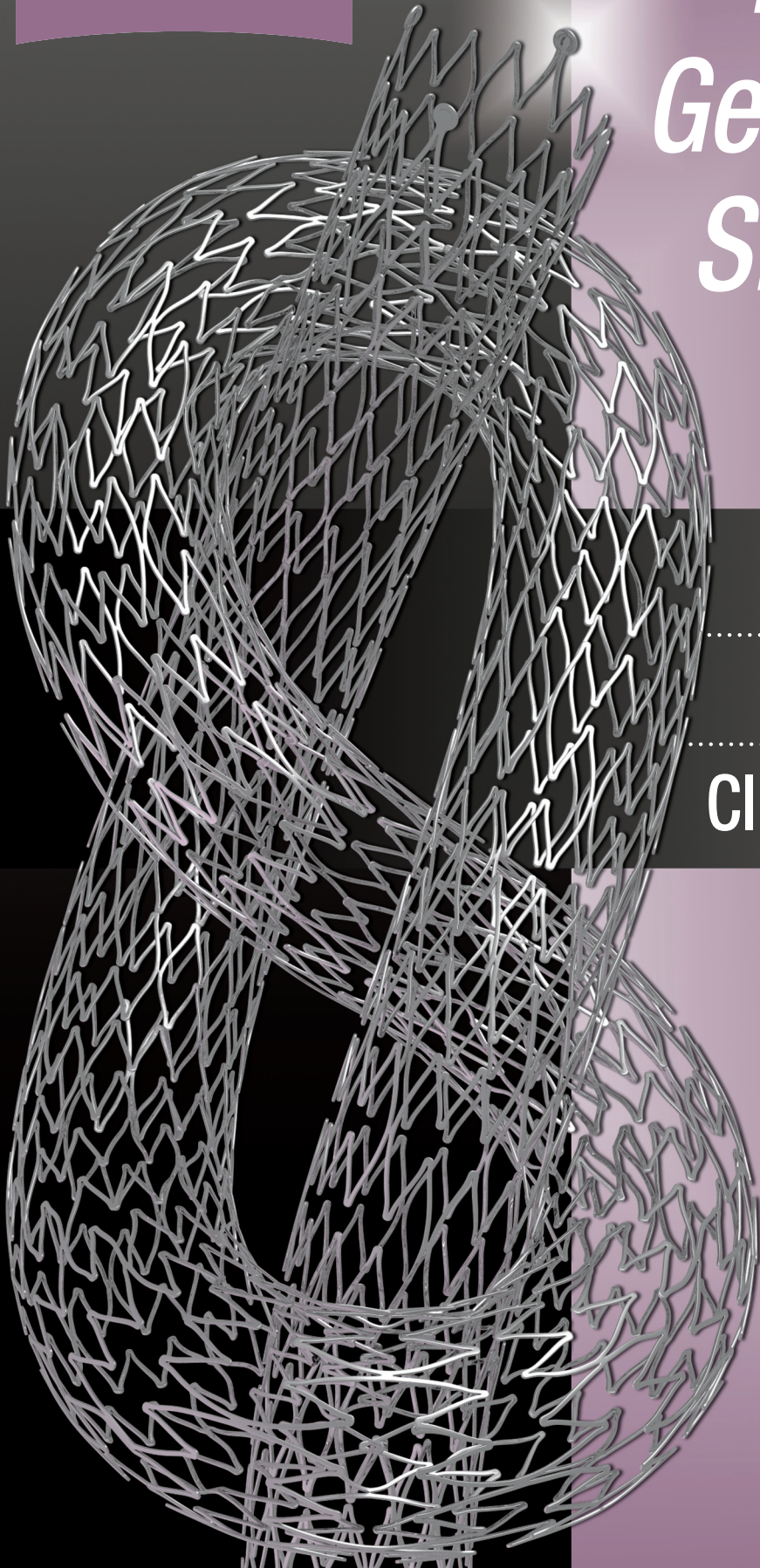

*The Next  
Generation  
SFA Stent*


Flexibility  
DURABILITY  
Clinical Evidence



True 2nd Generation stent design


EverFlex provides superb strength and excellent three-dimensional flexibility for challenging SFA stenting.

- Spiral connections significantly improve flexibility and vessel conformability without sacrificing radial force.
- Minimal shortening during deployment

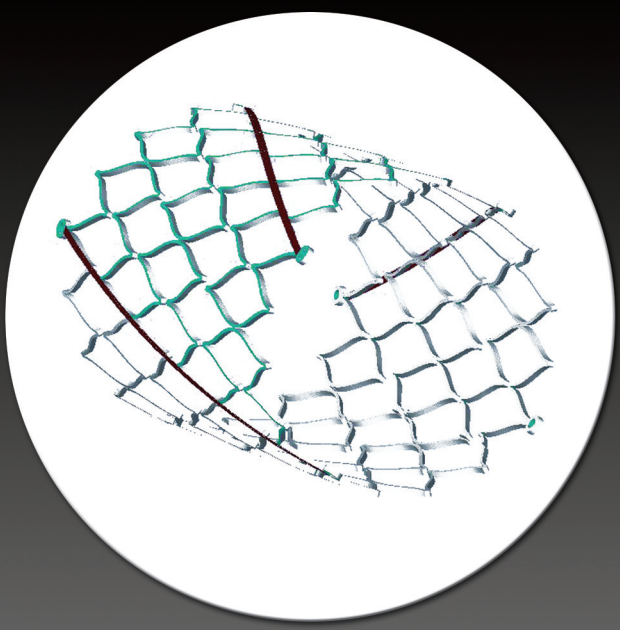

flexibility

durability

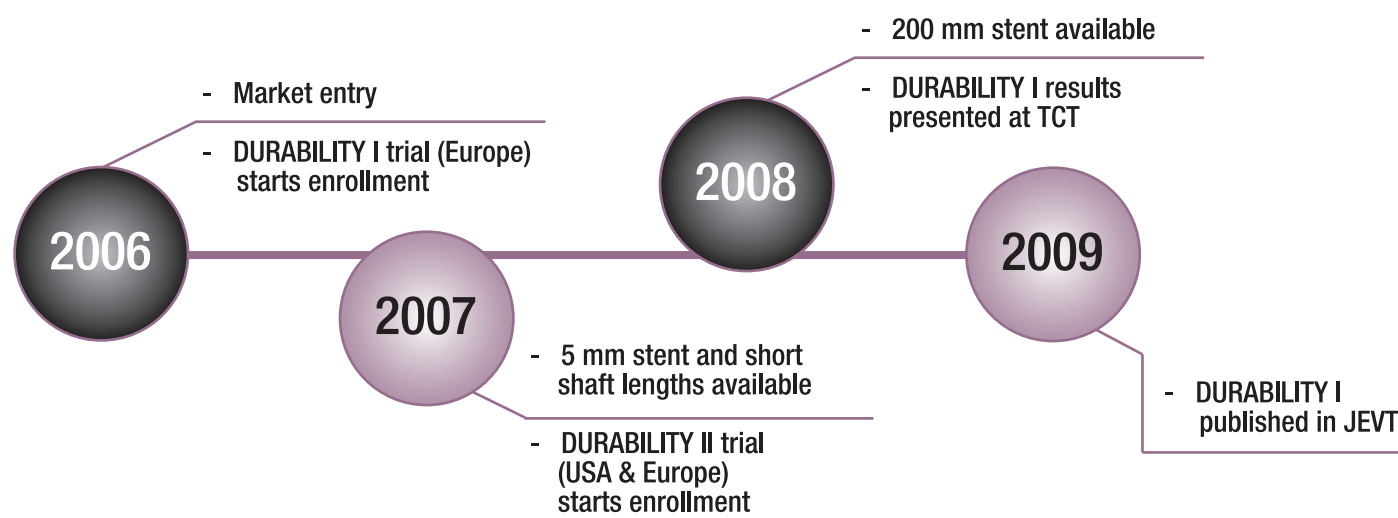

clinical evidence



Superior DURABILITY

More wave peaks between each interconnection distribute stress across multiple struts, greatly enhancing fatigue performance, durability and uniformity of deployment.

Stent manufacturing and finishing:

Superb fatigue performance through

- Laser Cutting
- Microblasting
- Electropolishing
- Surface passivation
- Riveted markers

**DURABILITY**

Study measuring the durability of the PROTÉGÉ EverFlex stent in lesions of the superficial femoral artery

Nitinol Stent Implantation in Long Superficial Femoral Artery Lesions: 12 MONTH RESULTS OF THE DURABILITY I STUDY  
Bosiers M, et al: Journal of Endovascular Therapy 2009;16(3):261-269

Results (n=151 patients)	
Lesion Length:	96.4 mm
Primary Patency:	72.2%
Freedom from TLR:	79.1%
Rutherford Class:	2.84 baseline to 0.60 at 1 year (significant)
ABI Improvement:	0.6 baseline to 0.9 at 1 year (significant)

141 out of 151 patients received only 1 single EverFlex stent.