




PRONTO[®] V4

Extraction Catheter

**Superior deliverability
and kink resistance**

Pronto V4 Features

- Embedded longitudinal wire enhances deliverability and kink resistance
- Uniformly large extraction lumen maximizes thrombus extraction
- Lubricious hydrophilic coating promotes vessel navigation and trackability
- Patented Silva distal tip provides vessel protection
- Full range of 5.5F, 6F, 7F and 8F sizes

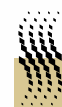

vascular
SOLUTIONS


Superior Deliverability

In simulated tortuosity testing, Pronto V4 provided maximum

DELIVERABILITY TESTING



Pronto V4

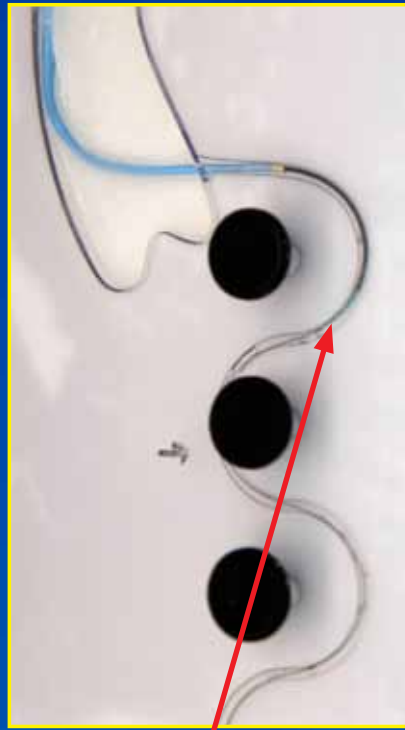


Manufacturer "E"




Manufacturer "F"

Uniformly Large Extraction Lumen

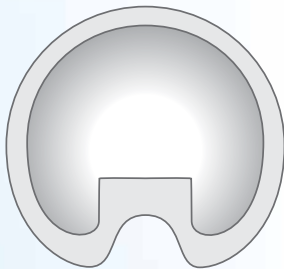

Proximal Section

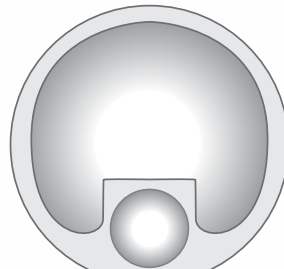

Distal Section

- Large full-length lumen provides a powerful aspiration force with no restrictions

Lubricious Hydrophilic Coating


- Distal 18cm hydrophilic coating promotes smooth navigation through tortuous anatomy
- Hydrophilic coating enhances trackability and reduces friction by maintaining lubricity throughout clinical procedure



y and Kink Resistance

imum distal deliverability with superior kink resistance.*

KINK RESISTANCE TESTING

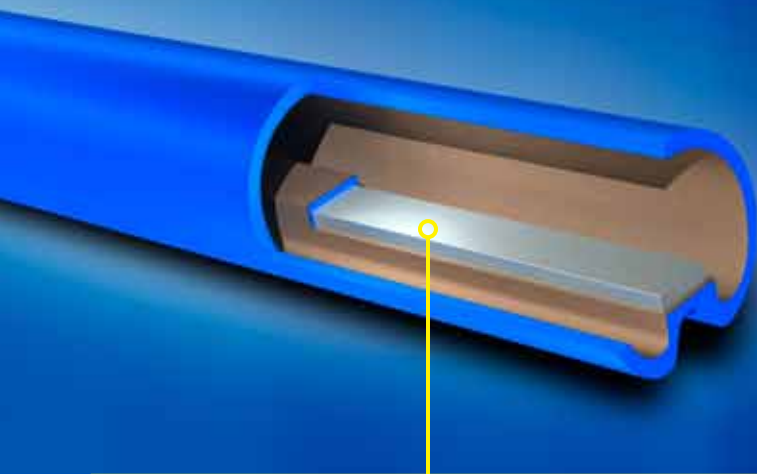


Embedded longitudinal wire in the Pronto V4 enhances deliverability and kink resistance



Pronto V4



Manufacturer "E"

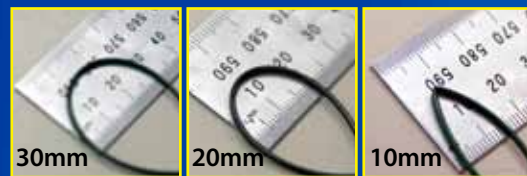


Manufacturer "F"

* Data on file


Patented Tapered Silva* Tip

- Tapered entry profile—designed to improve navigation and protect vessel wall during entry and aspiration
- Self-centering— designed to extract thrombus while preventing vessel wall adhesion
- Sloped extraction lumen opening—designed to aspirate larger thrombi

Pronto Silva tip

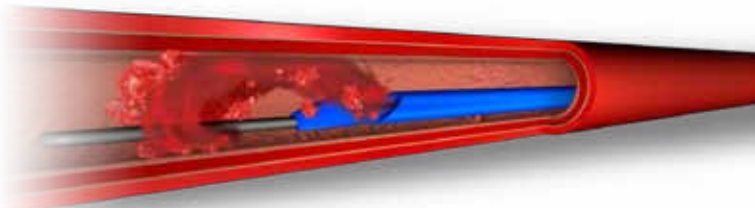

vs. Standard open-end extraction catheter tip

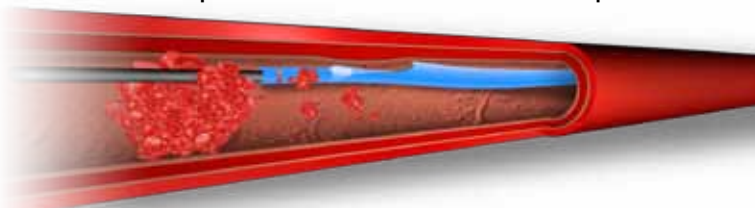

* Designed by Dr. Pedro Silva, Niguarda Hospital, Milan, Italy. U.S. and International patents.

PRONTO V4

Extraction Catheter

The Pronto V4 extraction catheter is indicated for the removal of fresh, soft emboli and thrombi from vessels in the coronary and peripheral vasculature.





Ordering Information*

Version	Model	Guide Catheter Compatibility	Sheath Compatibility	Guidewire Compatibility	Minimum Vessel Size	Working Length
5.5F	4005	6F ($\geq 0.066"$)	5F	0.014" (0.036mm)	1.75mm (0.068")	138cm
6F	4006	6F ($\geq 0.070"$)	5F	0.014" (0.036mm)	2.00mm (0.078")	138cm
7F	4007	7F ($\geq 0.078"$)	6F	0.014" (0.036mm)	2.25mm (0.088")	138cm
8F	4008	8F ($\geq 0.088"$)	7F	0.014" (0.036mm)	2.50mm (0.098")	138cm

*Packaged in quantities of 1 unit per box

Specifications:

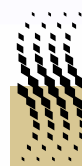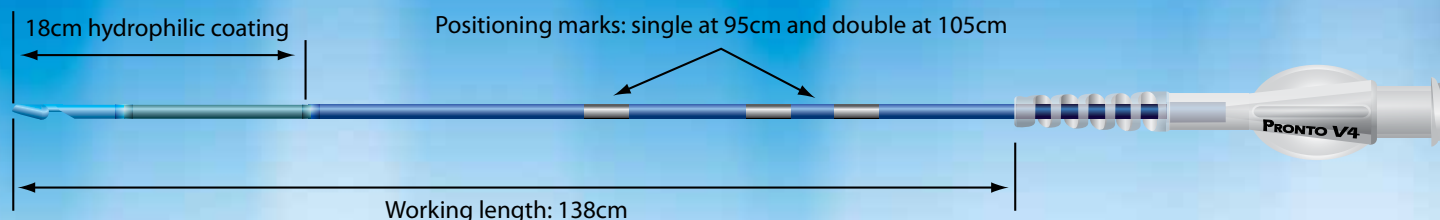

Vascular
SOLUTIONS

Vascular Solutions, Inc.
6464 Sycamore Court
Minneapolis, Minnesota 55369 USA

Customer Service:
United States: 888.240.6001
International: (001) 763.656.4298
customerservice@vascularsolutions.com

www.vascularsolutions.com

Please see the Pronto Instructions for Use for a complete listing of the indications, contraindications, warnings and precautions.

U.S.A. CAUTION: Federal Law (U.S.A.) restricts this device to sale by or on the order of a physician.

U.S. and International patents.

Pronto is a registered trademark of Vascular Solutions, Inc.

©2010 Vascular Solutions, Inc. All rights reserved. ML2315 Rev. A 12/10