

PRONTO[®]-SHORT

Extraction Catheter

Dialysis Graft and Native Fistula Declotting System



Make *SHORT* Work of Thrombus ... *Pronto!*

The Pronto-Short extraction catheter is indicated for:

- the removal/aspiration of embolic material (thrombus/debris) from vessels of the arterial system.
- the removal/aspiration of thrombus from clotted synthetic dialysis grafts and arteriovenous fistulas.

The Pronto-Short catheter is contraindicated in:

- vessels <2.5mm in diameter
- the removal of fibrous, adherent or calcified material (e.g., chronic clot, atherosclerotic plaque)

*Designed for percutaneous
thrombus aspiration
from dialysis grafts and fistulas*





Simple, one-operator deployment

- 45cm OTW catheter
- 80cm PTFE coated platinum tipped guidewire included
- Two locking 30ml aspiration syringes
- No powered components




The Pronto-Short delivers a maximized extraction lumen through a 6F introducer sheath with .018" guidewire compatibility.

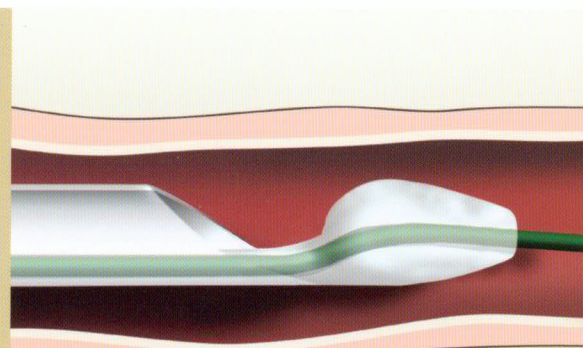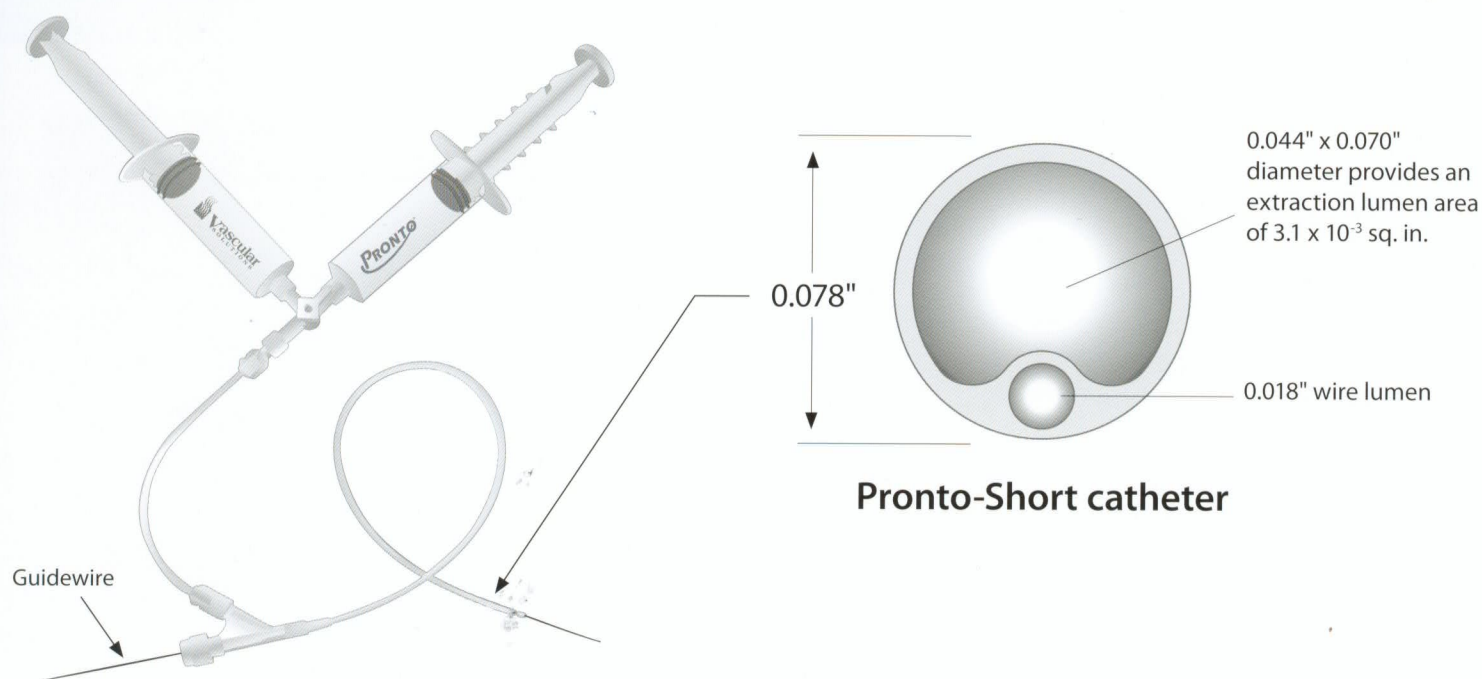

Proprietary Silva* tip for vessel protection

- Self-centering wire lumen to protect vessel wall during entry
- Designed to prevent extraction lumen contact with the vessel wall

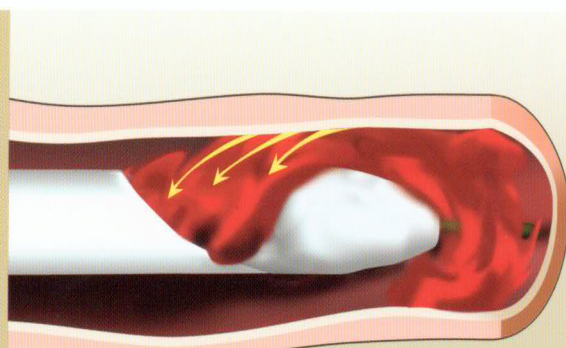

Measured extraction

- Graduated locking syringe
- 3-way stopcock and second locking syringe for added extraction capacity

PRONTO-SHORT

Extraction Catheter

Model	Description	Box Quantity
5030	Pronto-Short Extraction Catheter	1

Pronto-Short Specifications

Introducer sheath compatibility:	≥6F
Guidewire compatibility:	≤0.018"
Overall length:	44.5cm
Working length:	40cm
Marker band:	Located ≈3mm from the distal tip
Syringes:	Two locking 30ml polycarbonate syringes



SAN LUIS POTOSÍ No. 96 COL. ROMA SUR 06700 MÉXICO, D.F.
TEL.: 5574-9019 • FAX: 5564-9095
dimsamexico@mac.com
CALIDAD HUMANA®

U.S.A. CAUTION: Federal Law (U.S.A.) restricts this device to sale by or on the order of a physician.

CAUTION: The Pronto-Short Extraction Catheter procedure should be performed by physicians with adequate training in the use of the device.

Please see the Pronto-Short Instructions for Use for a complete listing of the indications, contraindications, warnings and precautions.

U.S. and International patents pending.
Pronto is a registered trademark of Vascular Solutions, Inc.

©2008 Vascular Solutions, Inc. All rights reserved.

ML1541 Rev. E 03/08



Vascular Solutions, Inc.
6464 Sycamore Court
Minneapolis, Minnesota 55369 USA
USA Customer Service: 888.240.6001
www.vascularsolutions.com

ECIREP

TheraGenesis GmbH
Bruehlstrasse 50
76297 Stutensee, Germany
Tel: +49 (0) 72 4995 2063
Fax: +49 (0) 72 4995 2064

CE
0344